



# SERVICE BULLETIN

## KODAK X-OMAT Processor

Eastman Kodak Company, Health Imaging, Rochester, New York 14650

SERVICE BULLETIN NO. 191


January 1994

**3133 KODAK M35A-M X-OMAT Processor**  
**3408 KODAK M35-M X-OMAT Processor**  
**3047 KODAK M35A X-OMAT Processor**  
**3041 KODAK M35 X-OMAT Processor**

### Installing the Receiving Bin on the M35 Series Processors

When installing the RECEIVING BIN, Part No. 261183, on the M35 series of processors, you may experience difficulties fitting the 2 BIN FASTENERS through the holes in the tabs and metal brackets. This is caused by the holes being misaligned.

The RECEIVING BIN can be exchanged using the Green Tag Procedure or the tip edge of the tab can be filed approximately  $\frac{1}{8}$  -  $\frac{1}{4}$  in. to obtain proper alignment of the fastener holes.

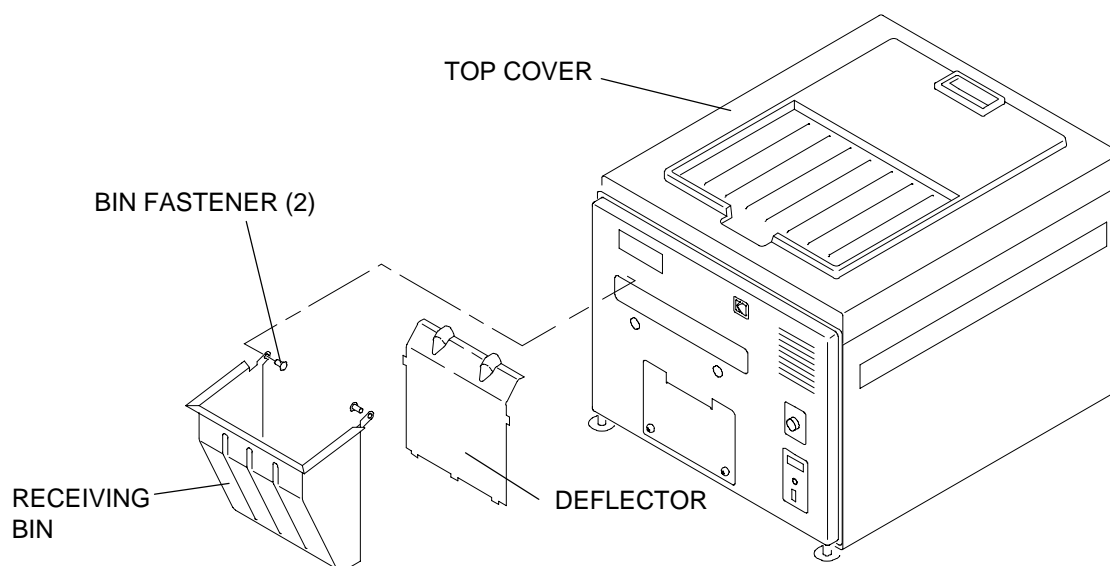

H112\_0151BCA  
H112\_0174BC

**For More Information:**

**U.S. Distributors Contact:**

Health Imaging  
Customer Support Operations  
1-800-336-4722

**CES Personnel Contact:**

Film Handling TAC  
U.S.: 1-800-822-1414  
Canada: 1-800-433-1414  
International: 716-724-4675

*Kodak*, and *X-Omat* are trademarks.

N-902

