




SITE SPECIFICATIONS
for the
Kodak X-Omat 460 RA Processor

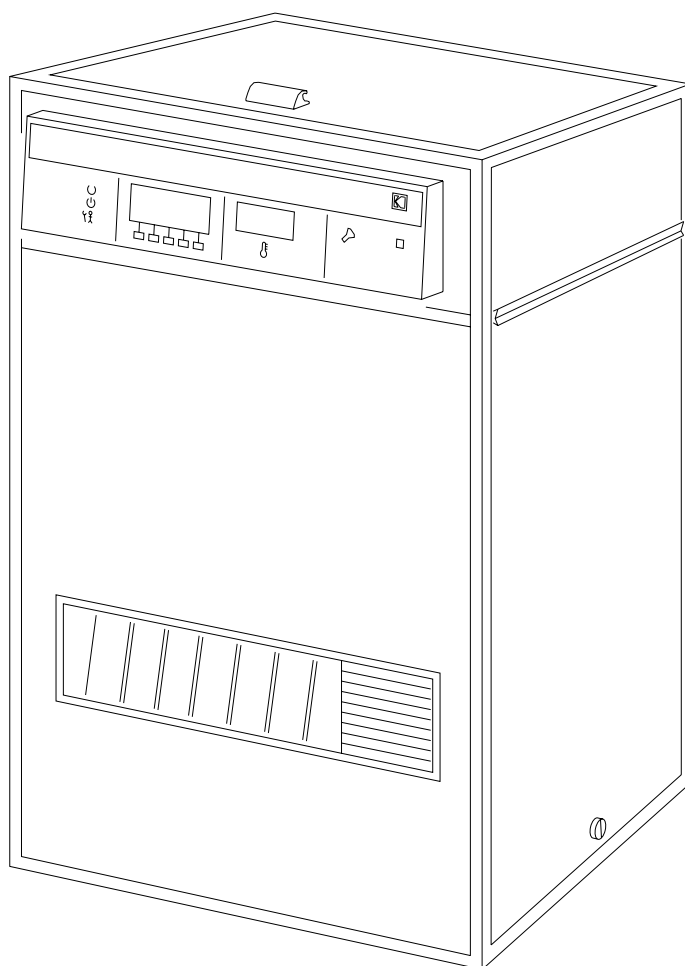

H108_0003DA

PLEASE NOTE

The information contained herein is based on the experience and knowledge relating to the subject matter gained by Eastman Kodak Company prior to publication.

No patent license is granted by this information.

Eastman Kodak Company reserves the right to change this information without notice, and makes no warranty, express or implied, with respect to this information. Kodak shall not be liable for any loss or damage, including consequential or special damages, resulting from the use of this information, even if loss or damage is caused by Kodak's negligence or other fault.

Table of Contents

Description	Page
Overview of the Processor	3
Dimensions	5
Electrical Requirements	10
Basic Requirements.....	10
Standard Service Options.....	10
Processor Interface Connector.....	12
Water and Drain Requirements	13
Water Supply	13
Drain.....	13
Main Power Disconnect (Wall mounted, not furnished).....	14
Environmental Requirements	16
Room Ambient Conditions.....	16
Air and Heat.....	16
Center of Gravity	18

Overview of the Processor


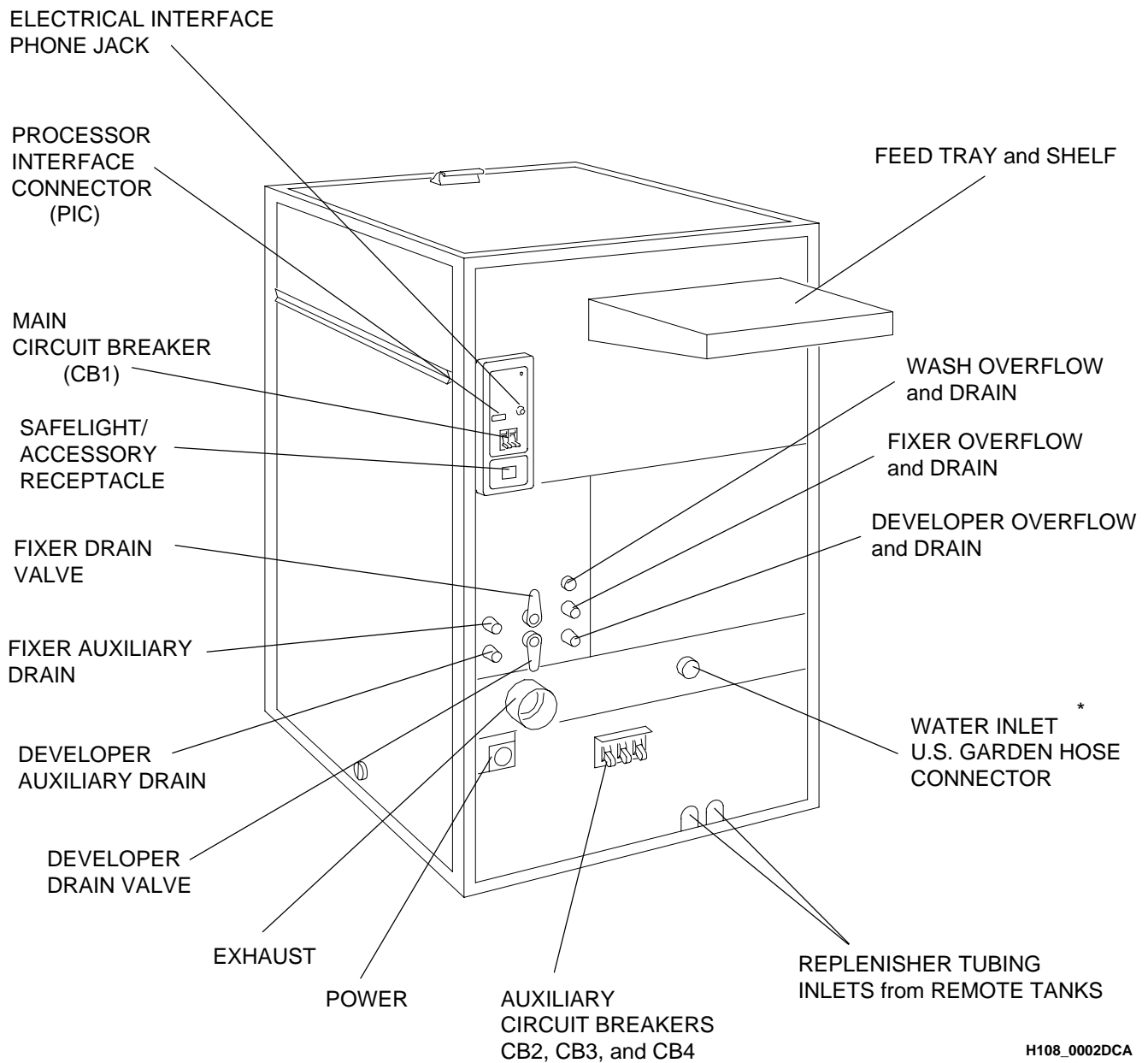



Figure 1 Feed-End View

* Supplied in the pre-pack is an adapter for 1/2-inch NPT.

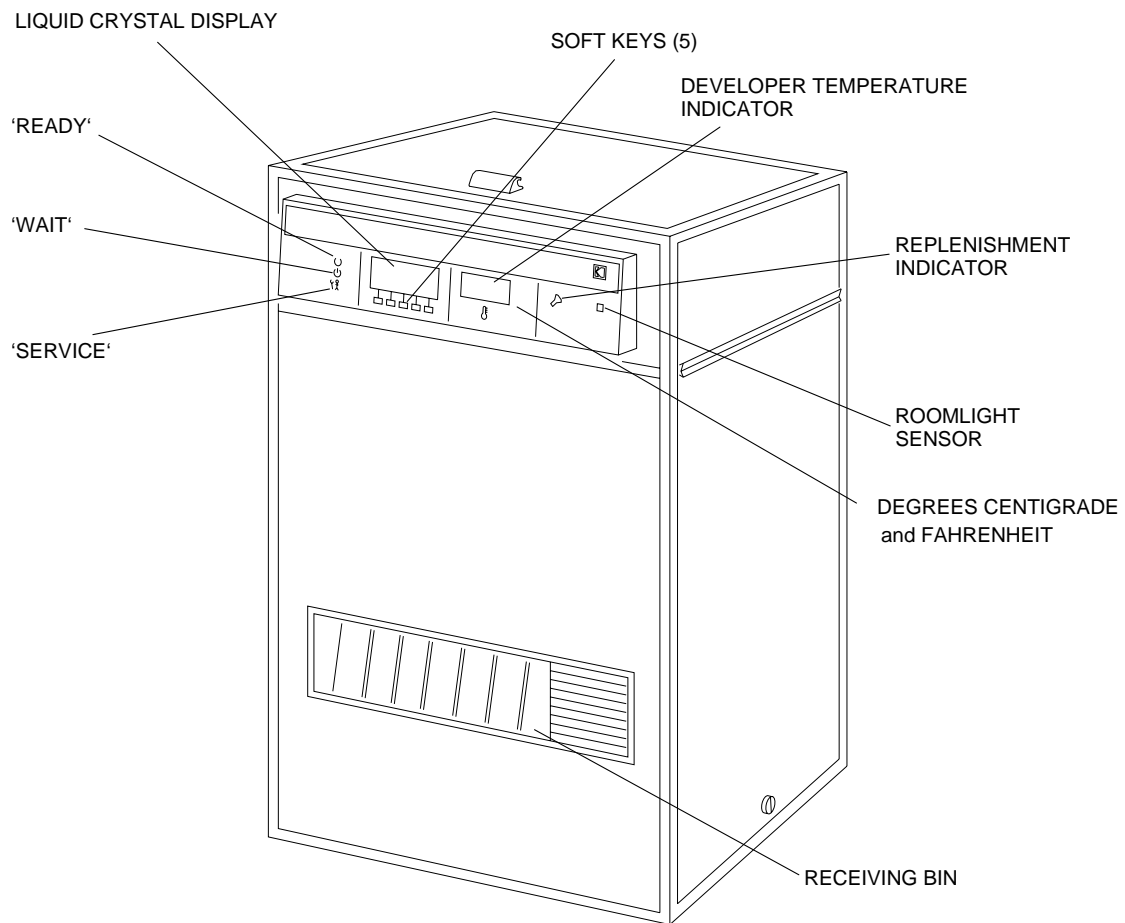

H108_0003DAA
H108_0003DA


Figure 2 Receiving-End View

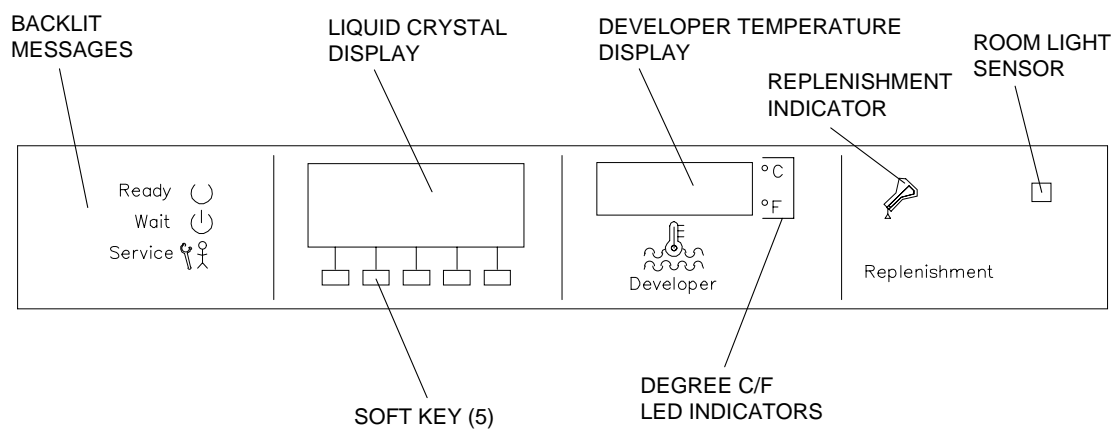

H104_0080BCA
H104_0080BA

Figure 3 Display Panel on Receiving-End of the Processor

Dimensions

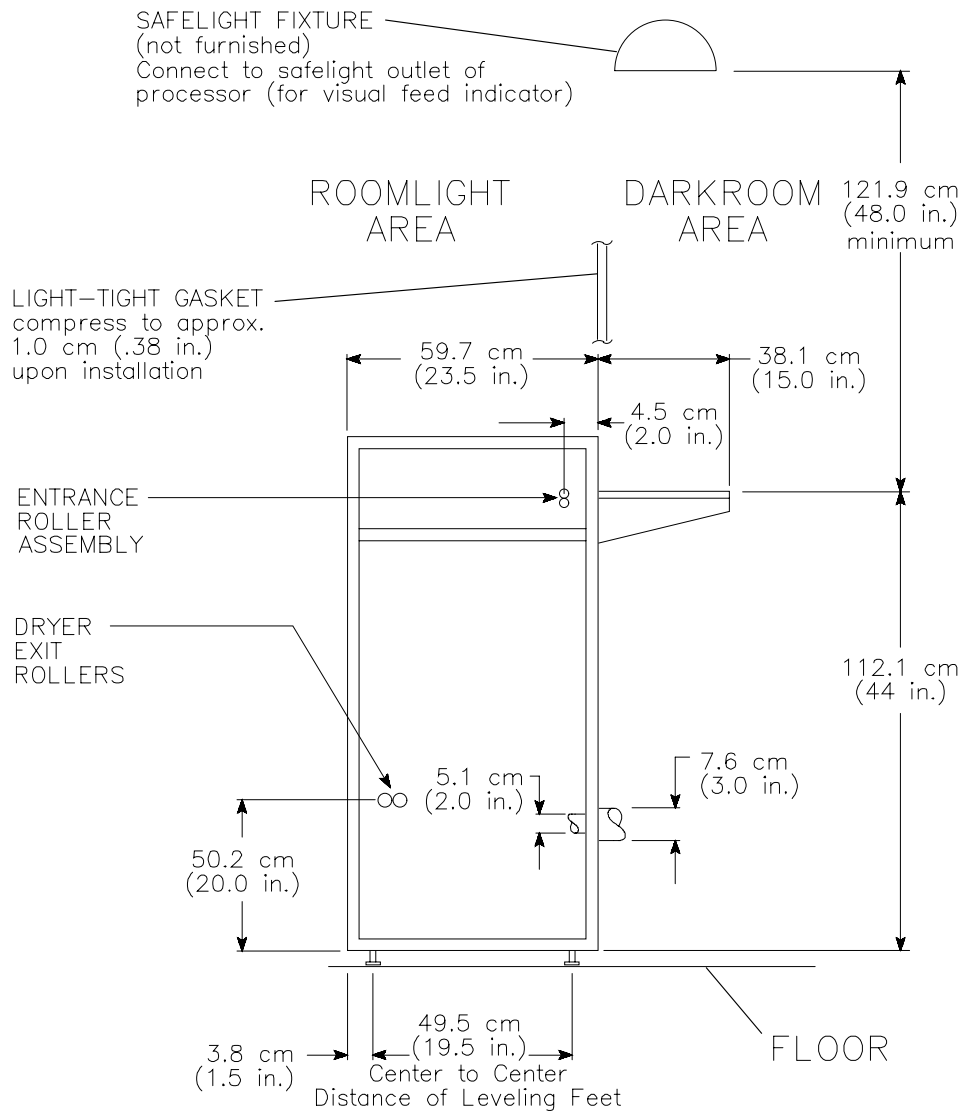


H108_0007DA

Figure 4 Side Dimensions

Table 1 Dimensions and Weight of the Processor

Description	Crated	Uncrated
Length	73.7 cm (29 in.)	66.0 cm (23½ in.) Without feed tray 97.8 cm (38½ in.) With feed tray
Width	88.9 cm (35 in.)	76.2 cm (30 in.)
Height	141.0 cm (55½ in.)	123.2 cm (48½ in.)
Weight (Tanks Empty)	214 kg (471 lb)	201 kg (422 lb)
Weight (Tanks Full)	Not Applicable	235 kg (499 lb)

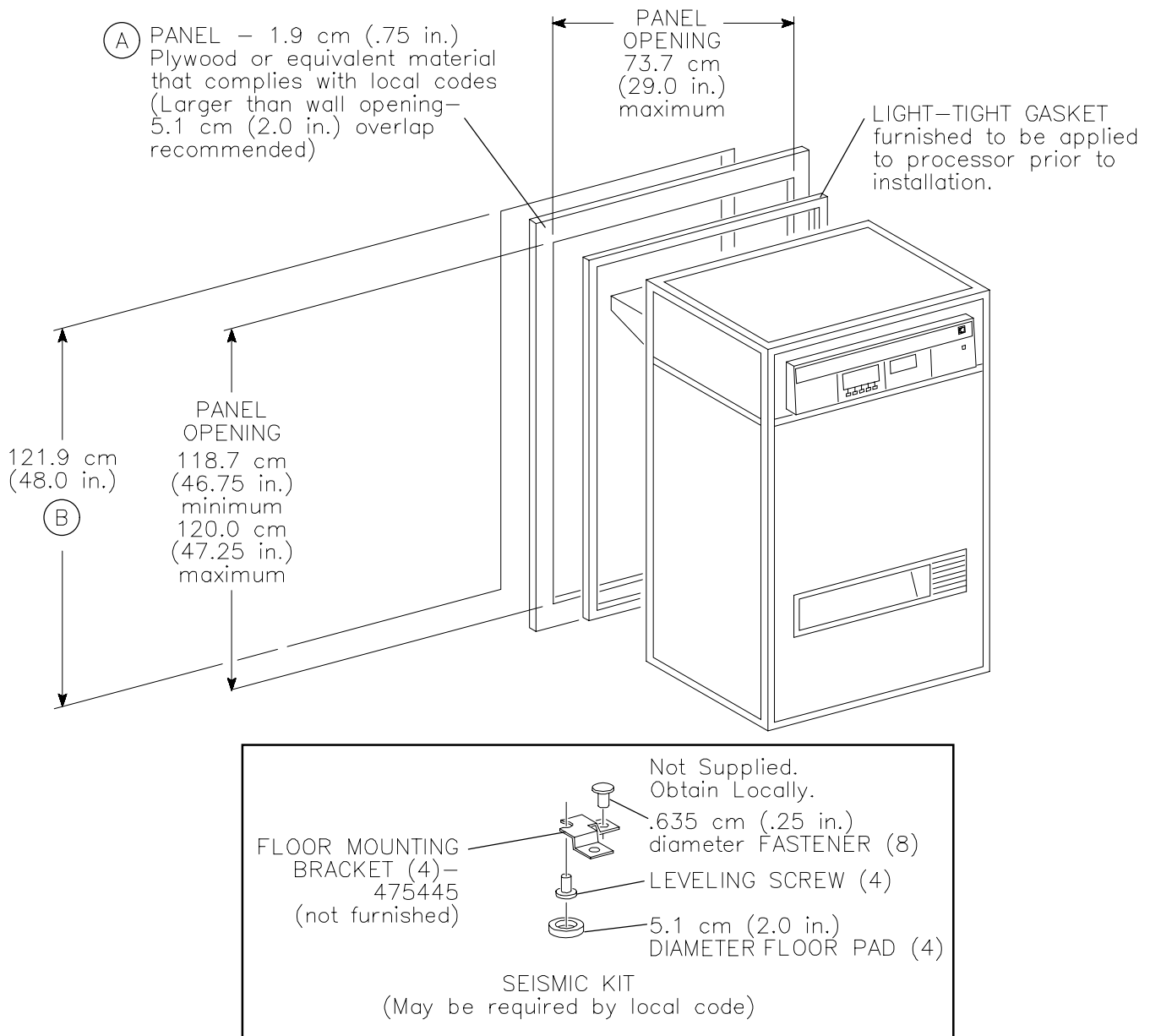


H108_0010DA

Figure 6 New Wall Installation - Feed End

- A** If the wall around the opening is straight and exactly perpendicular to the floor, this panel may not be necessary. The wall opening dimensions should match the inside dimensions of the panel opening.
- B** Make sure that the vertical dimension of 48±1/16 in. for the wall opening is measured from the finished floor.

Replenishment Tank Requirements

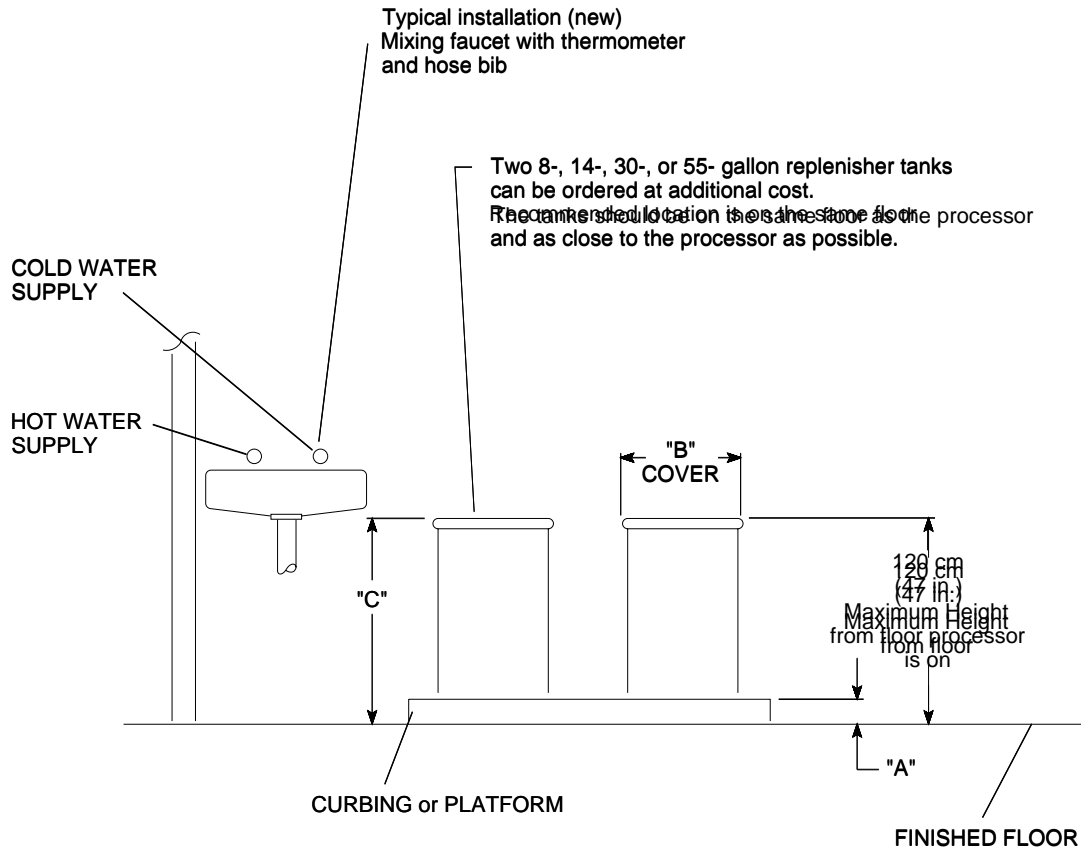


H108_0005DAA
H108_0005DA

Figure 7 Replenisher Tanks

Table 2 Dimensions for Replenisher Tanks

DESCRIPTION	DIMENSION	14 GAL	30 GAL	55 GAL
Max Platform Height	"A"	48.3 cm (19 in.)	35.6 cm (14 in.)	15.2 cm (6 in.)
Tank Diameter	"B"	43.2 cm (17 in.)	55.9 cm (22 in.)	61.0 cm (24 in.)
Tank Height	"C"	58.4 cm (23 in.)	70.5 cm (27¾ in.)	90.8 cm (35¾ in.)
External Replenishment Tank Area	"D" x "E" (MIN)	61.0 X 127.0 cm (24 X 50 in.)	61.0 x 152.4 cm (24 x 60 in.)	66.0 x 172.7 cm (26 x 68 in.)

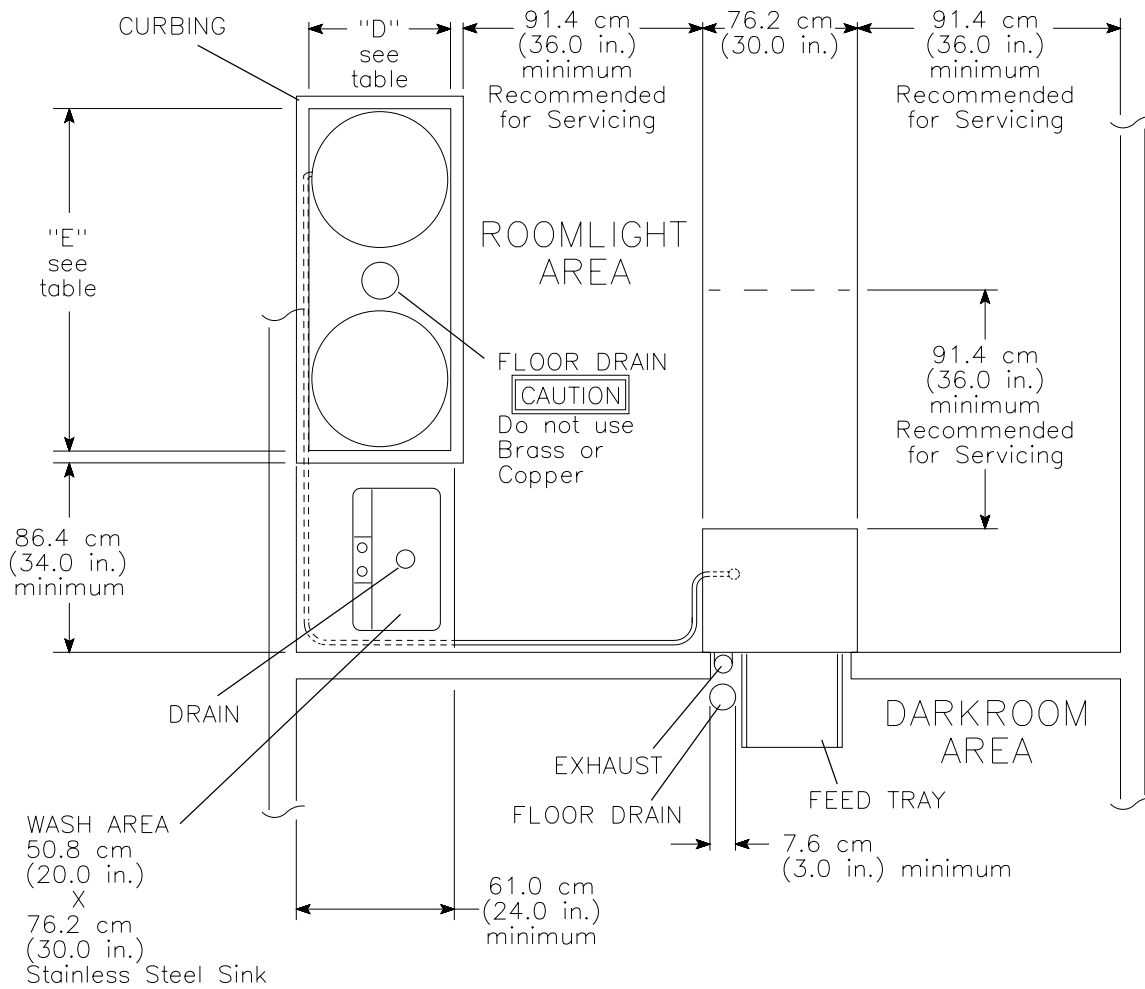

H108_0006DA

Figure 8 Suggested Room Layout Showing Replenishment Tanks and Drain Locations

Table 3 Maintenance and Operation Access Requirements

Description	Recommendation
Receiving End of Processor	91.4 cm (36 in.)
Feed End of Processor	91.4 cm (36 in.)
Drive Side of Processor	91.4 cm (36 in.)
Non-Drive Side of Processor	91.4 cm (36 in.)
Top of Processor	91.4 cm (36 in.)

Electrical Requirements

- Basic Requirements**
- 35 A, single-phase
or
 - 25 A, three-phase
 - single-phase 2-wire, or 3-phase 3 or 4-wire service
 - earth ground


IMPORTANT

All electrical services, *including* **earth ground**, must comply with local and national electrical codes.

Standard Service Options

Table 4 Service Options

Voltage Volts	Frequency Hz	Service
200	50/60	2 - Wire Single Phase plus Earth Ground
220	50/60	2 - Wire Single Phase plus Earth Ground
230	50/60	2 - Wire Single Phase plus Earth Ground
240	50/60	2 - Wire Single Phase plus Earth Ground
100/200	50/60	3 - Wire Single Phase plus Earth Ground
120/240	50/60	3 - Wire Single Phase plus Earth Ground
200	50/60	3 - Wire 3 - Phase Delta plus Earth Ground
120/208	60	4 - Wire 3 - Phase Wye plus Earth Ground
127/220	50	4 - Wire 3 - Phase Wye plus Earth Ground
220/380	50	4 - Wire 3 - Phase Wye plus Earth Ground
230/400	50	4 - Wire 3 - Phase Wye plus Earth Ground
240/415	50	4 - Wire 3 - Phase Wye plus Earth Ground

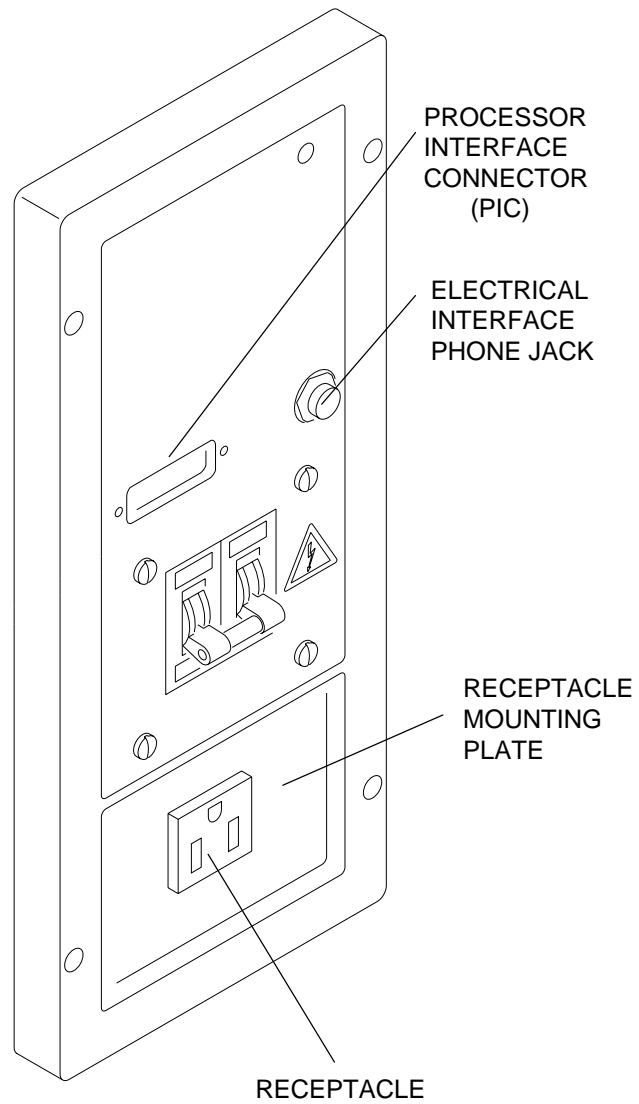

H108_0141CCA
P108_0141C

Figure 9 Control Panel

NOTE

The RECEPTACLE MOUNTING PLATE is removable for use of other receptacles.

Processor Interface Connector PIC Connector on Processor: 15 Position Panel Mount Subminiature D, AMP Part No. 1-747299-0.

Customer Supplied Mating Connector: AMP Part No. 205206-1 or Equivalent.

Table 5 Pinout for the Processor Interface Connector (PIC) Port

PIC#	Signal Description	Type	Function
1.	Shield	Ground	Connection for the cable's shield in an RS232 system.
2.	Transmit Data (TxD)	Output	Data out from the processor in an RS232 system.
3.	Receive Data (RxD)	Input	Data into the processor in an RS232 system.
4.	Request to Send (RTS)	Output	Handshake line for the control of data out from the processor.
5.	Clear to Send (CTS)	Input	Handshake line for the control of data into the processor.
6.	Common/Return	Input	Signal common in an RS232 system and +5v return.
7.	Processor Cycle C0	Output	TTL logic signal indicating processor cycle ¹ . 0 is ≤ 0.5 V @ 4mA 1 is ≥ 3.84 V @ 4mA
8.	Processor Cycle C1	Output	TTL logic signal indicating processor cycle ¹ . 0 is ≤ 0.5 V @ 4mA 1 is ≥ 3.84 V @ 4mA
9.	$\overline{\text{Run}}$	Input	TTL logic signal, active LOW input for placing processor in Run mode.
10.	Alarm	Output	TTL logic signal, active HIGH indicating an alarm condition.
11.	Service	Output	TTL logic signal, active HIGH indicating a fault requiring service.
12.	Ready	Output	TTL logic signal, active HIGH indicating processor is operating within specifications.
13.	$\overline{\text{Film Feed}}$	Output	TTL logic signal, active LOW indicating a sheet of film may be inserted into the processor.
14.	Reserved for Future Use	Input	TTL logic signal input not used at this time.
15.	+5 Volt DC (1 A Max.)	Output	Logic supply source of 5 V dc also used to indicate processor is on.

¹ Cycle	C1	C0
K/RA	1	1
Rapid	1	0
Standard	0	1
Extended	0	0

Interface Phone Jack Specification

Accepts standard 3 circuit 1/4 in. phone plug:
SWITCHCRAFT No. 290 or equivalent

Water and Drain Requirements

Water Supply

- Temperature: 4°C (40°F) to 32°C (90°F)
- Pressure: 173 to 448 kPa (25 to 65 psi); Install regulator if required.
- Flow: Controlled within the processor to 5.7 L/min (1.5 gal/min), $\pm 10\%$
- Filtration: 50 micron filter required

IMPORTANT

Water supply must comply with local codes; **do not** use iron piping.

NOTE

- If the upper limit of the room ambient temperature or the water supply temperature is exceeded, the developer and fixer temperatures may not be controlled correctly.
- Tempered water service is suggested for cleaning the processor and for mixing chemicals manually.

Drain

- Capacity: 15 L/min (4 gal/min)
- Connection: Open drain; avoid a solid connection to the processor to eliminate suction conditions


IMPORTANT

Do not use brass or copper for the drain lines. Drain service must comply with local codes.

Main Power Disconnect (Wall mounted, not furnished)

A main power disconnect switch, as shown in the Figure, consisting of a 2-pole for single-phase and 3-pole for 3-phase, thermomagnetic circuit breaker with solid neutral and common trip **or** a fused disconnect switch **must** be—

- (1) located on a wall adjacent to the processor in the lighted area.
- (2) easily accessible from the processor site.
- (3) visible from the processor site.

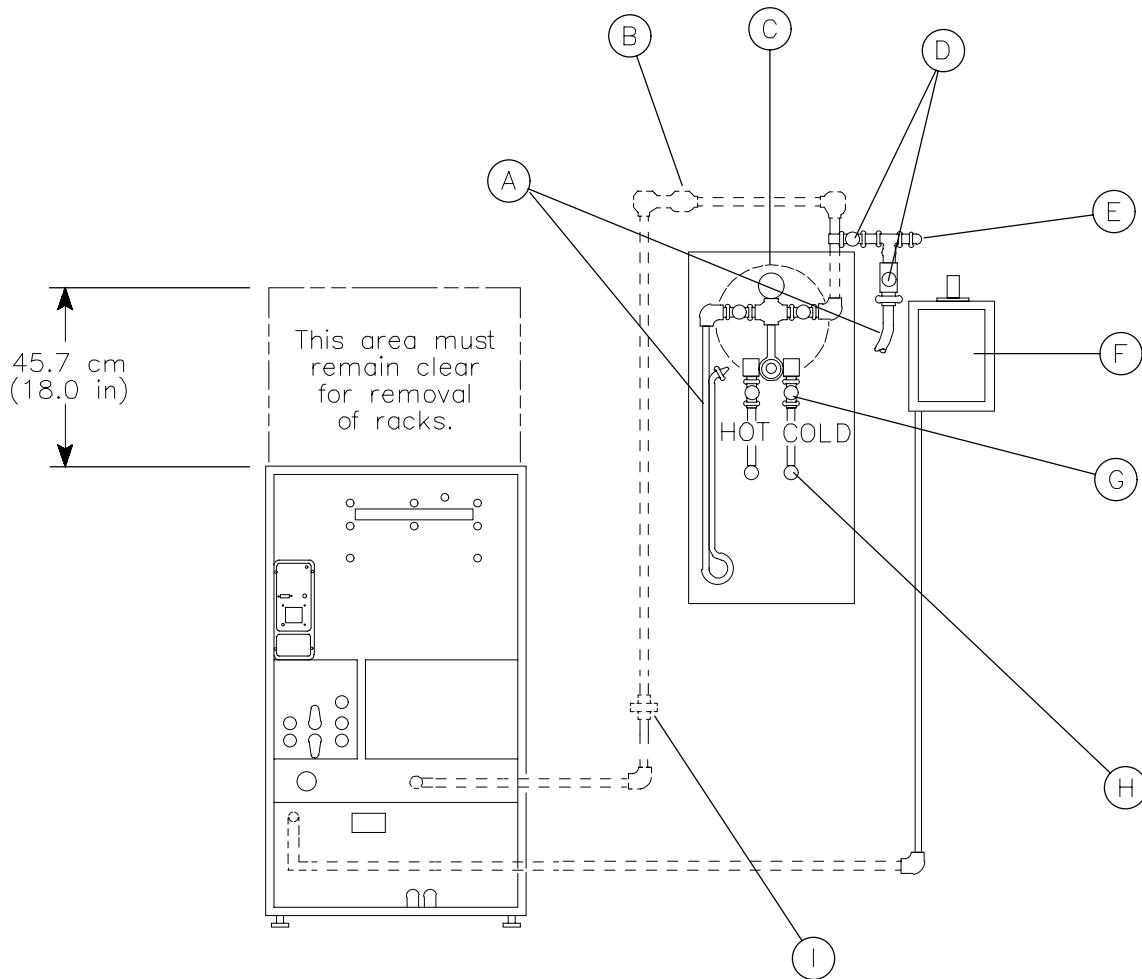

H108_0009DAA
H108_0009DA

Figure 10 Electrical and Water Connections

Table 6 Electrical and Water Connections (cont.)

A	Service Hose - suggested length to reach the Processor and Replenisher Tanks.
B	13 mm (½ in.) NPT Check Valve
C	Kodak Thermostatic Mixing Valve 13 mm (½ in.) NPT available from Eastman Kodak Company, Part No. 467621 (Not supplied with the Processor.)
D	Shutoff Valves 13 mm (½ in.) NPT. Two additional required. Available from Eastman Kodak Company, Part No. 459981 (Not supplied with the Processor).
E	Cold Water Supply 13 mm (½ in.) NPT
F	35 AMP - 2 POLE for single phase, or 25 AMP - 3 POLE for 3-phase service, Thermomagnetic Circuit Breaker. (Locate safe distance from water service.) (Not supplied with the Processor.)
G	Shutoff Valves 13 mm (½ in.) NPT. Two additional required. Available from Eastman Kodak Company, Part No. 459981.
H	Hot and Cold Water Supply 13 mm (½ in.) NPT.
I	13 mm (½ in.) NPT Union. Locate as close to the Processor as possible.

IMPORTANT

A mixing valve is not required for the processor if incoming water temperature is between 4° and 32°C (40° and 90°F).

Follow local electrical and plumbing codes.

NOTE

The processor is UL listed, CSA certified, and TUV approved. Specifications are subject to change without notice.

Pass service through the wall to the feed end of the processor in darkroom. Service controls may be located on either side of the processor for easy accessibility.

Environmental Requirements

Room Ambient Conditions

- Temperature: 16°C (60°F) to 30°C (86°F)
- Humidity: 15% to 76% relative humidity, non-condensing

Air and Heat

Air Exhaust (full load):

- 75 CFM maximum
- Temperature: 66°C (150°F) maximum
- Moisture: 600 grains/min

Heat Load to Room:


- 2 Kw (6824 BTU/hr or 1719 kcalories)

Building Exhaust Duct:

- Use AIR METER TL-2431 and modified J TUBE (CHECK TUBE 592380) to measure negative static pressure in the EXHAUST DUCT 30.5 cm (12 in.) from the end that is to be connected to the processor. See Figures 11 and 12.

NOTE

The processor must be off and the duct must not be connected while you are making the measurements.

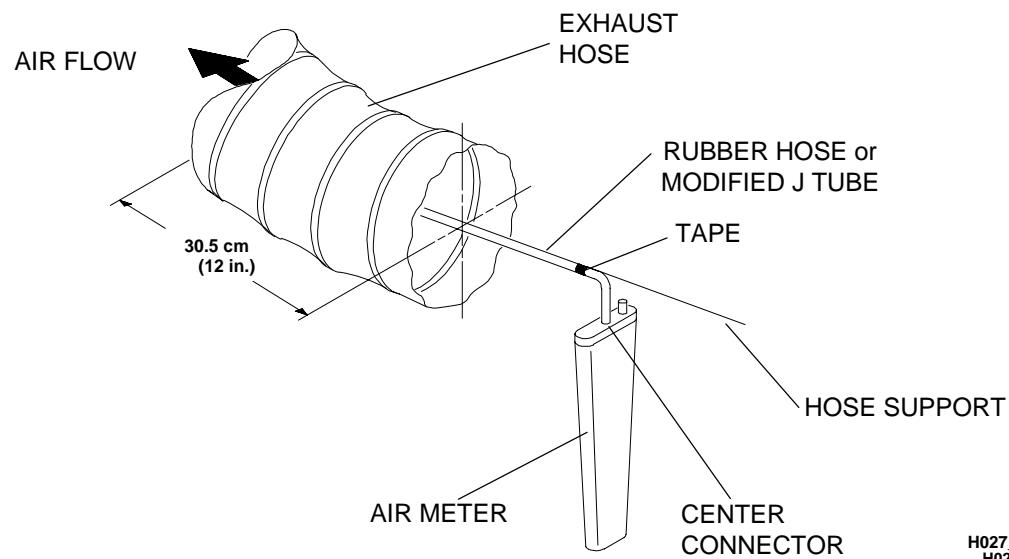


H027_0100BCA
H027_0100BA

Figure 11 Measuring the Static Pressure

- Adjust the clearance between the building exhaust duct and the duct from the processor as shown in Figure 12 to obtain the required static pressure outlined in Table 7.

Table 7 Required Static Pressures


Duct Diameter	Negative Static Pressure, (Water Head)	
	MIN	MAX
7.6 cm (3 in.)	0.76 mm (0.03 in.)	1.02 mm (0.04 in.)
10.2 cm (4 in.)	0.25 mm (0.01 in.)	0.51 mm (0.02 in.)

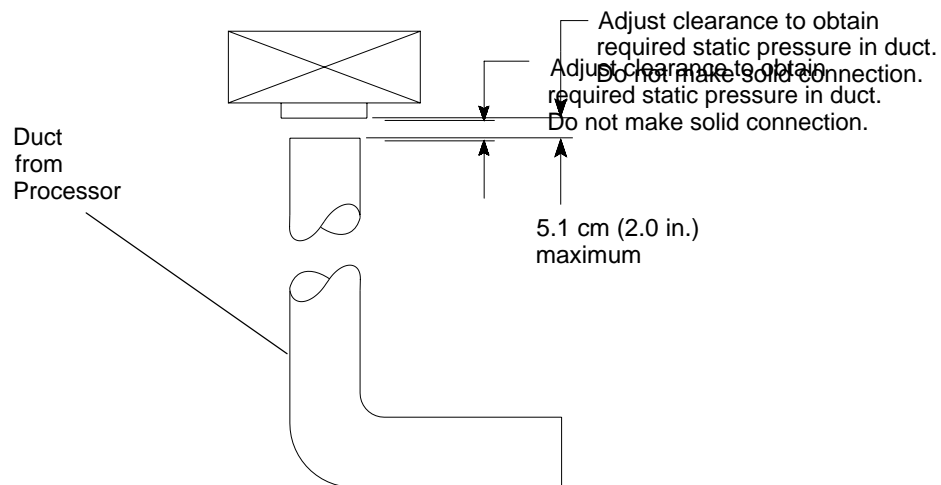

H108_0008BAA
H108_0008BA

Figure 12 Exhaust Requirements

Center of Gravity

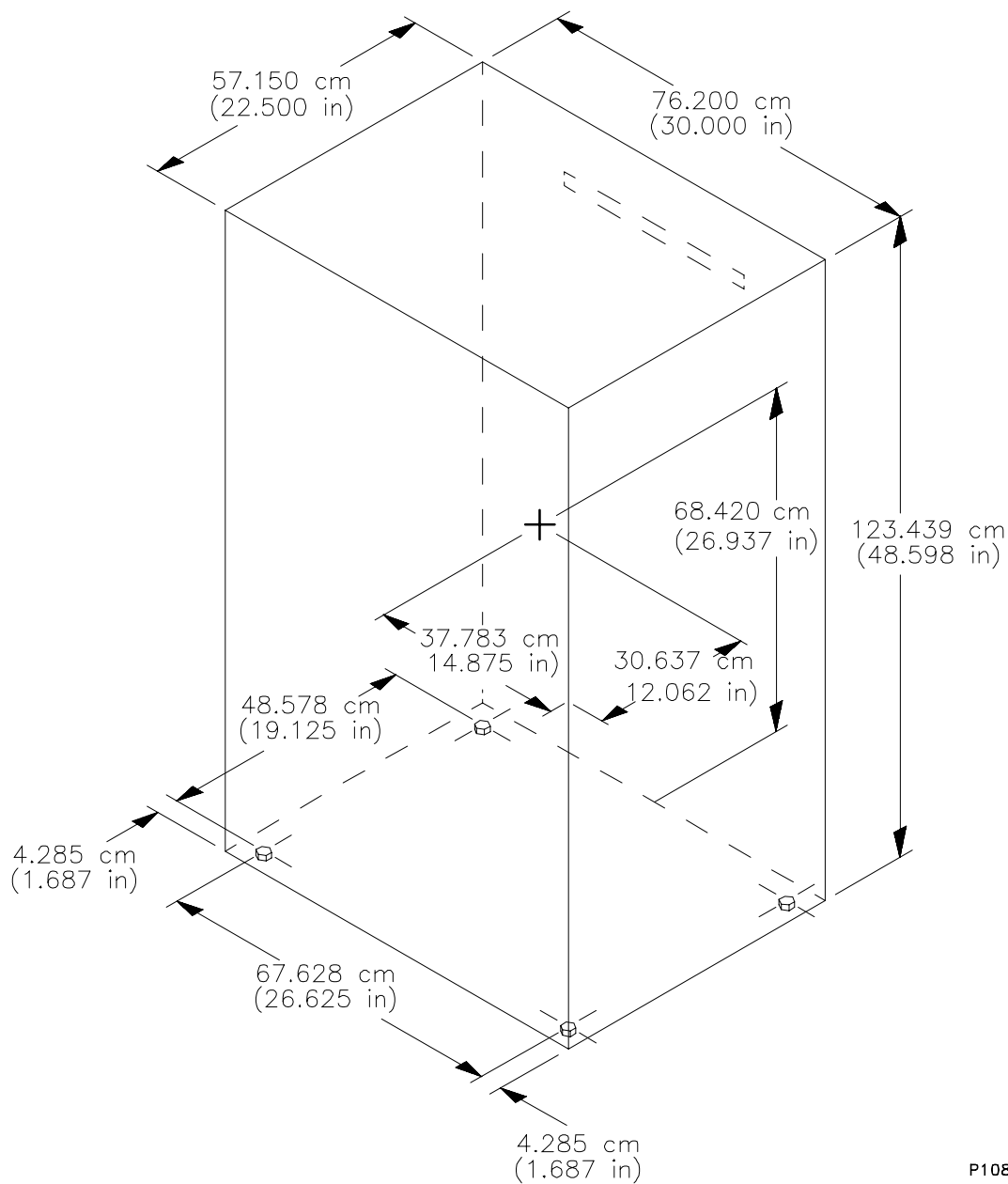

P108_0086DA

Figure 13 Center of Gravity

NOTE

Center of gravity is shown for the processor with the DEVELOPER, FIXER, and WASH TANKS empty.

BLANK PAGE

3056ss_c.txt

Kodak and X-Omat are trademarks.

Printed in USA

EASTMAN KODAK COMPANY • ROCHESTER, N.Y. 14650

Health Sciences Division

